


ACTUATION SYSTEMS for NURSING HOME PATIENT LIFT APPLICATIONS

TiMOTION provides a simpler actuation solution for patient lifts. We offer an alternative for medical environments that require safe and easy-to-use patient transfer.

This solution ensures safety and comfortability in patient lifts at nursing homes.

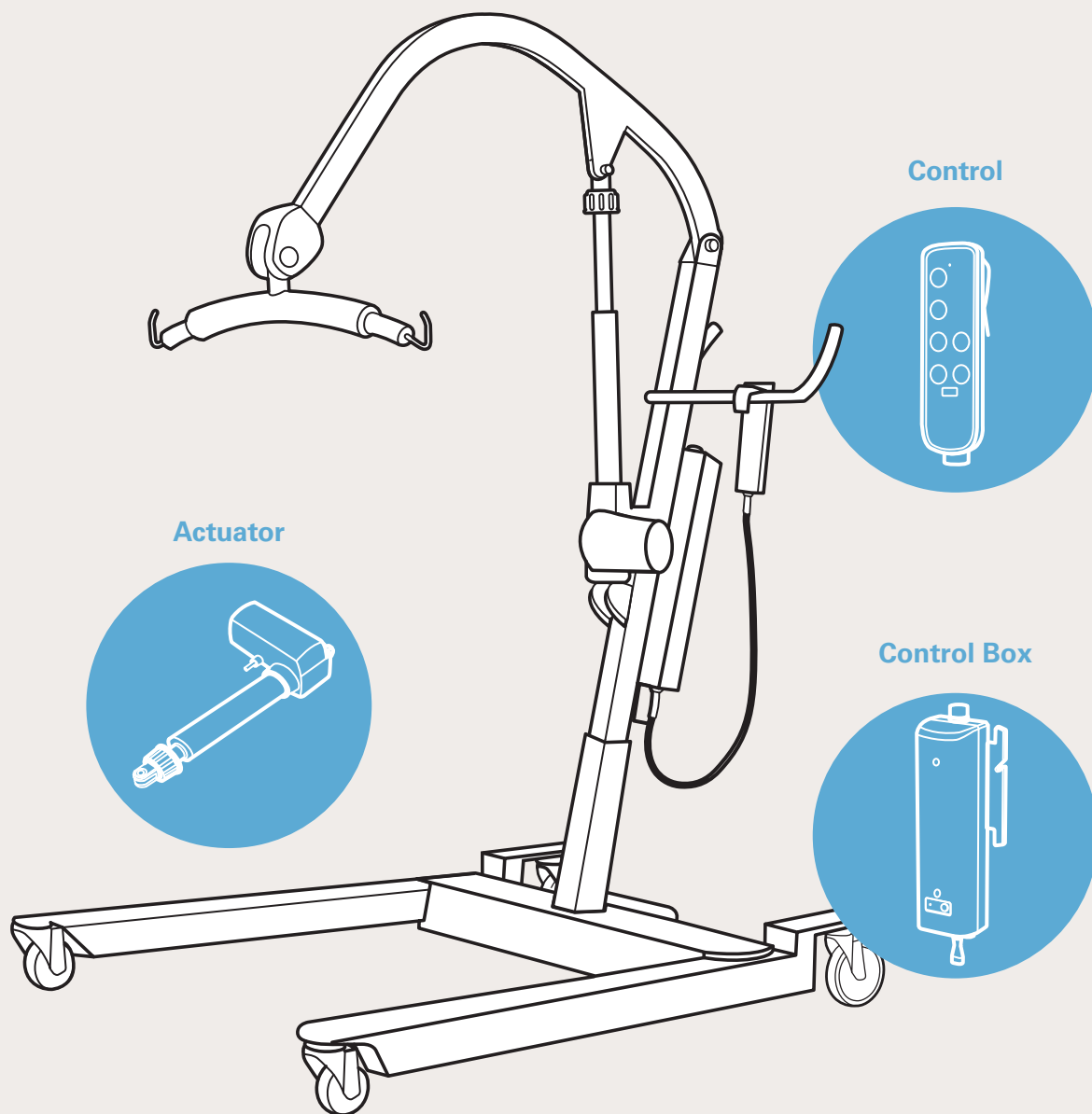

A SMART AND EASY TO USE ACTUATOR SYSTEM FOR PATIENT LIFTS

TC20 is a simpler control box for a patient hoist system which meets specific patient needs and secures a safe and comfortable transfer. It is portable and easily installable within nursing home applications.

TC20 KEY FEATURES:

- Built-in battery pack: 24V, 2.9Ah or 5Ah
- Charging with AC cable or plug directly
- IP rating: Up to IP66
- Cable management bracket (optional)
- Emergency stop button
- Easy installation, only two screws required
- LED power indicator on the battery pack (optional)


System Overview

Control box: TC20


Control: TH10, TMH3 or TMH13

Actuators for lifting: TA23, TA36 or TA37

Actuators for leg adjustment/ sling angle: TA16 or TA29

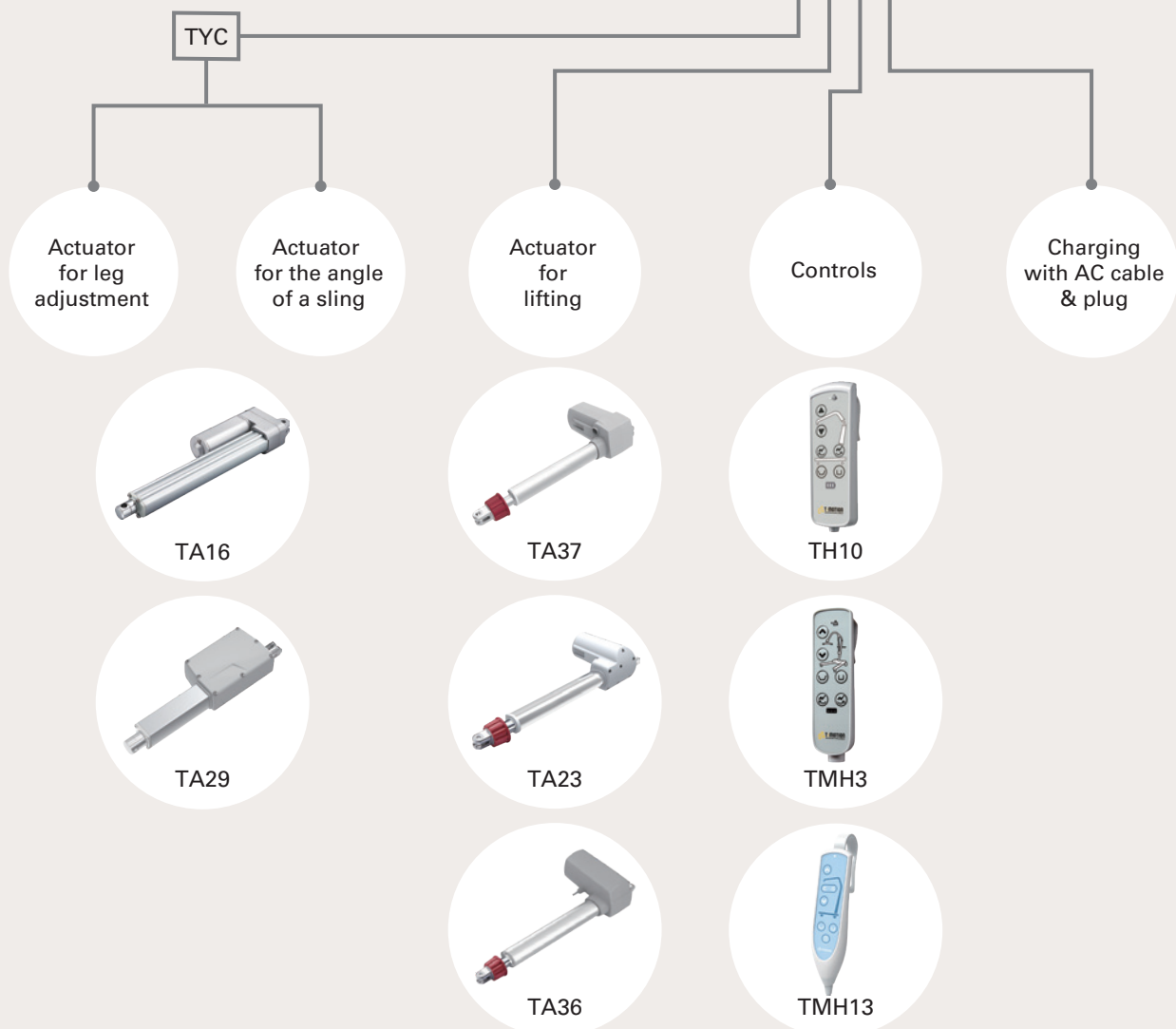

Actuator for Lifting

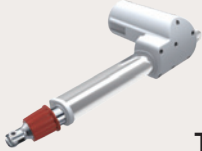Model	Maximum Load in Push	Maximum Speed at Full Load	IP Rating	Manual Release
 TA36	10,000N	5.4mm/s	IP66W	Available
 TA23	10,000N	3.2mm/s	IP66W	
 TA37	12,000N	5.3mm/s	IP66W	

Actuator for Leg Adjustment /Angle of Sling

Model	Maximum Load in Push	Retracted Length
 TA29	6,000N	≥ 178mm (depending on chosen options)
 TA16	3,500N	≥ Stroke+112mm



TMH13



TMH3




TH10

CONTROL TMH13, TMH3 or TH10 SERIES


- TMH3 has a detachable handset cord
- Supports a maximum of three actuators
- LED indicator for the battery status
- Interface for lifting/ lowering, leg adjustments, and overload warning

Installation Guide

1. Install TC20 on the hoist frame
2. Install fixed screws of TC20. Make sure the screws are tightened properly

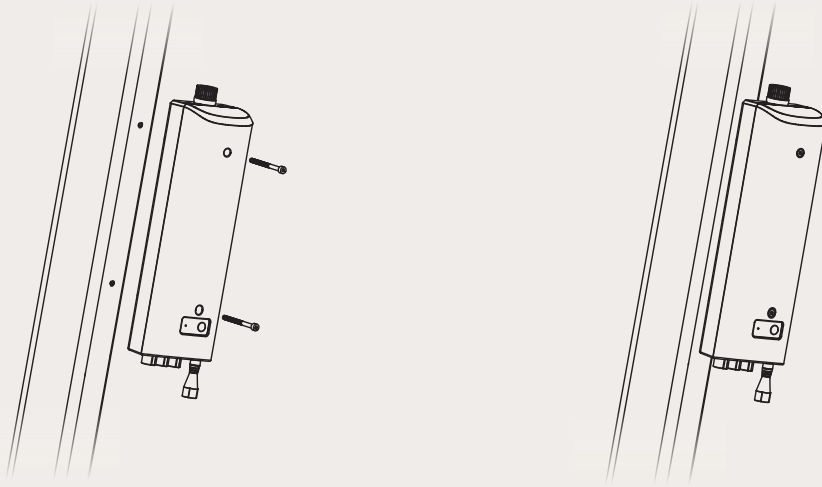

Standard Dimensions (mm)

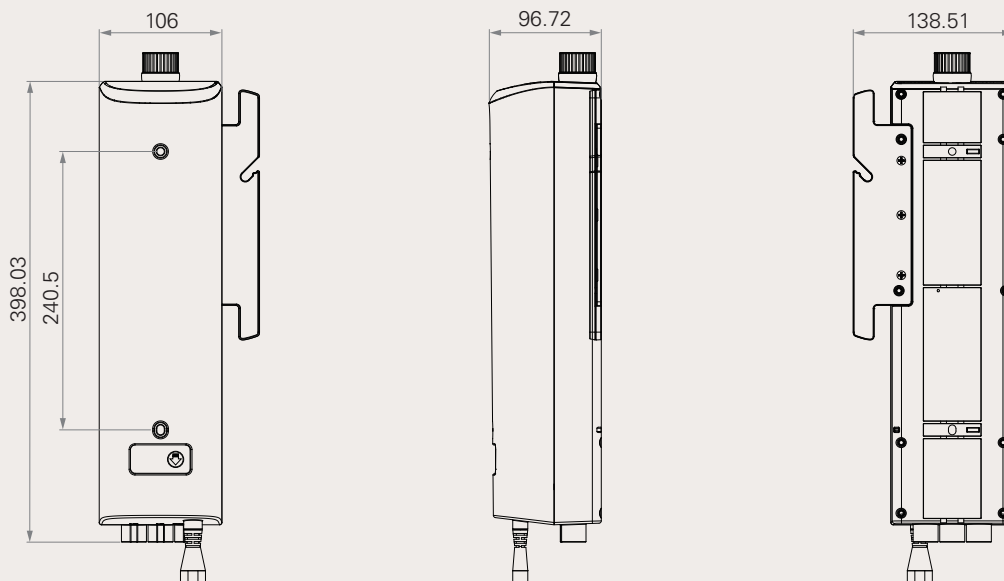

Taiwan & China

+886 2 2219 6633
sales.tw@timotion.com

Europe

+33 (0)1 74 82 50 51
sales.eu@timotion.com

USA

(704) 708 6924
sales.us@timotion.com

Australia

+61 3 8725 0331
sales.tw@timotion.com

Japan

+81 6 6713 1188
sales.jp@timotion.com

Korea

+82 31 745 1060
sales@timotion.co.kr

Russia

+7(495) 185-58-40
info@stalfors.ru